



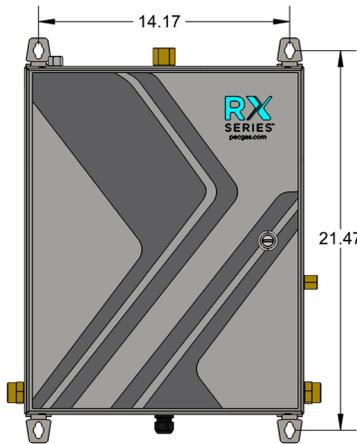

The **RX Series**® uses modern technology and industry leading UL / cUL certification to provide the medical gas market with a comprehensive control system that brings focus to customers' needs. The X-Treme Manifold Control® provides real time alerts to any electronic device, and updates on gas status and historical usage online, on an easy-to-read dashboard. This dashboard is configured to be viewed on a mobile device, offering flexibility when using the system from a phone, tablet, or PC. Our complimentary service is new to the industry, and truly a comprehensive alternative to other telemetry options.

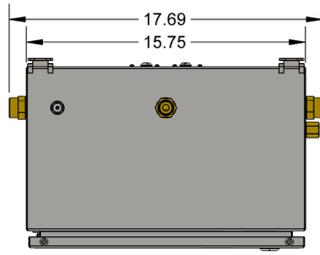

FEATURES

24/7 LEAK DETECTION

Will continuously monitor the reserve supply and provide leak indication lights if pressure is lost.

CYLINDER MANAGEMENT

Select from high-pressure or low-pressure source, including three different relief valve settings.

INTERNATIONAL POWER INPUT

AC 85 – 264V 50/60Hz 0.7A Class II (Includes a 6' Power Cord with NEMA 5-15P plug).

FULL-AUTOMATIC CHANGEOVER & EASY RESET

Once in-use side is depleted, the reserve side will flow uninterrupted, and the system will automatically switch to the reserve supply maintaining the integrity of the delivery gas supply. Manifold will reset the empty bank automatically on cylinder changeout.

GAS PRIORITY CONTROL

Can select a side being the primary supply when using a low-pressure and high-pressure mixed cylinder configuration by prioritizing the left or right bank.

3000 PSI

MAXIMUM INLET

AUTOMATIC ECONOMIZER

When low-pressure source is selected, the system will monitor and capture excess pressure on the reserve supply before reaching the programmed relief valve setting for each source.

DELAYED SWITCHOVER

In the event the in-use bank reaches the switchover point with low-pressure cylinders, the reserve side will flow uninterrupted for 20 minutes, before declaring an empty alert.

WIRELESS & LAN CONNECTIVITY

Connect to web-based telemetry using Wi-Fi network.

SERIAL 232/485 CONNECTION

Connect to Modbus or stream data using other protocols.

FLEXIBILITY

Uses individual low and high-pressure cylinders, and multi-cylinder packs interchangeably without modification.

WEATHERPROOF

Rated for IP65 location and can be placed outside in direct water spray or dusty conditions.

GAS STATUS LED'S

3-LED status indicators for each gas bank.

1- Economizer "ON" LED for each gas bank

1- "Leak Detected" LED for each gas bank

ACCURATE PRESSURE DISPLAY

Improved inlet pressure reading and switchover setting to 1-PSI increments.

SECURE CONTROLLER

Tamperproof user settings guarded by password protects the integrity of the controller from unauthorized use.

APPLICATIONS

Continuous Flow Without Power

Provide a continuous flow of gas to patients and equipment where uninterrupted flow is critical.

During a power loss, dual normally open solenoid valves will continue flowing gas from both primary inlets.

Hospital • Surgery Centers • Acute Care

Perfect for small or large health centers who adhere to the NFPA 99 source of standards.

Delivery Pressure

Models up to 200 psig delivery depending on gas type.

Hyperbaric Chambers

Supply pressurized oxygen to patients for decompression sickness and serious infections in either a room or a tube for treatment.

Nursing Facilities • Hospice

Providing a steady and controlled flow of oxygen for patient use within the facility.

Veterinary Hospitals

Equine Hospitals

MATERIALS

Mechanical Enclosure

Steel powder coated enclosure with latch, Rated IP65

Controller Enclosure

Polycarbonate with lockable latches, Rated IP65

Regulators

Brass bar stock and fabric reinforced neoprene diaphragms encapsulated PTFE seat

Tubing

3/8" stainless steel tubing with 2-ferrule tube connections.

Solenoid Valves

Brass 24 VDC coils reduce heat exchange and damage to seals.

Shutoff valves

Forged brass and conform to CGA V-9 and CGA V-1

Headers

Brass, hydrostatically tested to Minimum 4500 psig.

Pigtails

Flexible Stainless-Steel Teflon® or Tefzel® inner core, or rigid copper Cleaned, Capped, & Bagged

Cleaning

Meets CGA G-4.1

SYSTEM INCLUDES

UL cUL Listed:

RX Series® Gas Panel
X-Treme Manifold Control®

Cylinder Connections

Flexible stainless steel / rigid cooper pigtails with gas specific CGA connections on both ends. CGA specific check valve adapters are included at header inlet.

Pipeline Outlet Assembly

Includes pipeline relief safety relief valve set to either 200 or 500 psig depending on the final delivery pressure.

Headers

Expandable heavy wall brass. Rated 3000 PSIG working pressure

Shutoff Valves

Header master valves provide a layer of protection when replacing high-pressure cylinders



Mounting Hardware

Wall mounting brackets for gas panel, controller, and headers. Floor stand models include Mounting posts, post bases, header brackets with hardware

Cylinder Restraints and Chain

Cylinder restraints and chain will also be included for individual high-pressure gas sources.

Unistrut channel conforms to ASTM A1011 SS GR 33 E-Coat finish

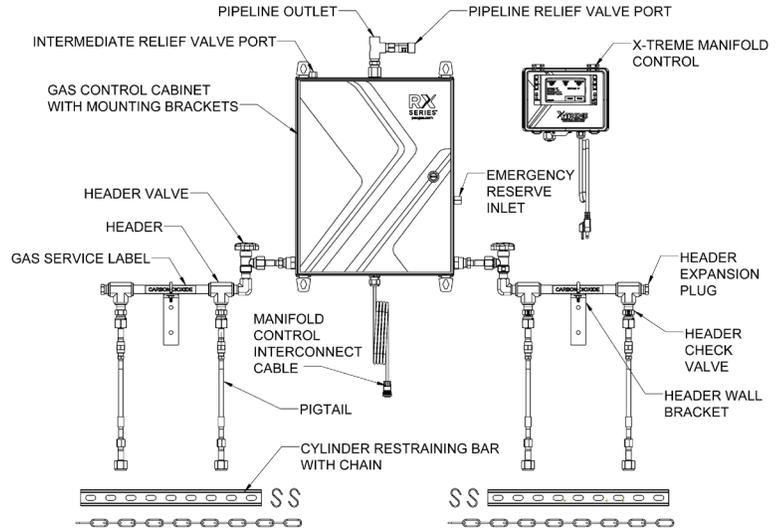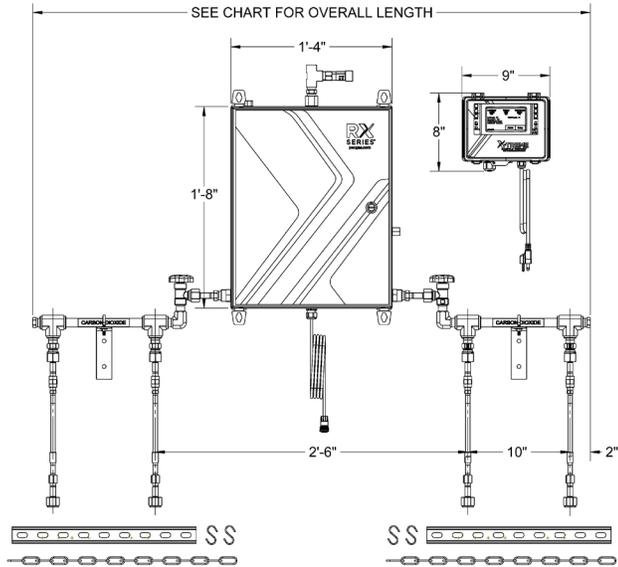

UL / cUL SERIES	GAS SERVICE	DELIVERY RANGE	RATED INLET PRESSURE
RX-1	Oxygen, CO2, Nitrous Oxide, Breathing Air	5-125 PSIG	3000 PSIG
RX-2	Nitrogen	5-200 PSIG	

NUMBER OF CYLINDERS & OVERALL LENGTH						
Header Type	2	4	6	8	10	12
Standard	2' - 4"	4' - 6"	6' - 2"	7' - 10"	9' - 0"	10' - 8"
Close Space	2' - 0"	4' - 0"	4' - 0"	4' - 10"	5' - 8"	6' - 6"